
J BIMBI

Slim

Manuale d'uso
User manual
Bedienungsanleitung
Manuel d'utilisation

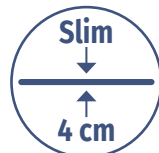12383


IT

EN

DE

FR

Display LCD / LCD display / LCD-Anzeige

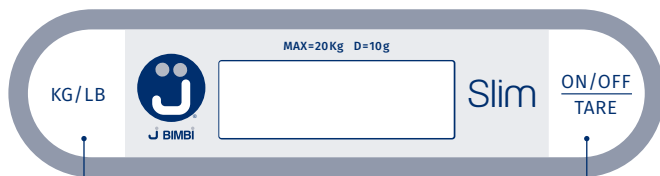Tasto "Unità di peso"
"Weight unit" button
Taste "Gewichtseinheit"
Touche "Unité de poids"

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Digital scale

Thank you for choosing a J BIMBI® product. This Slim digital scale is a useful and reliable tool for monitoring the growth of children right from the first days of life because it detects even the slightest changes in weight. The ergonomic and ultra-slim design makes it practical and easy to handle, while the automatic weighing lock function and the large backlit LCD display make it easy to read and detect the weight.

BEFORE USING THIS DEVICE READ THESE INSTRUCTIONS CAREFULLY AND STORE THEM FOR FUTURE REFERENCE.

RECOMMENDATIONS FOR USE

Place the scale on a flat, stable surface; uneven or soft surfaces (e.g. rugs, beds, sofas) may cause inaccurate weight measurements.

Make sure the surface of the scale is clean and dry before measuring your baby's weight.

HOW TO USE THE DEVICE

How to install the batteries

This scale operates with 4 standard AA batteries (included). Turn the scale upside down and install the batteries in the battery compartment making sure to respect the polarity as shown. Close the cover and place the scale on a flat, stable and solid surface.

When the batteries are low the display shows «**LO**» (1), open the battery compartment and replace the batteries. Remove the batteries if you plan to leave the scale unused for a prolonged period of time.

kg/lb weight unit setting button

Located to the left of the display, this button allows you to choose the unit of measure when weighing: kilograms (kg) or pounds (lb). Select the desired weight unit when the scale is on and «**0.00**» is displayed (2) (otherwise the unit of measure cannot be changed).

How to use the scale

1. To obtain accurate measurements place the scale on a flat, stable surface. Uneven or soft surfaces can cause inaccurate weight measurements.
2. Press the "ON/OFF - TARE" button to power on the scale. The display will show «**Hello**» (3) and after a few seconds «**0.00**» (2).
3. Press the "Weight unit" button to select the desired weight unit (kg or lb).
4. Gently place the baby on the scale, paying special attention to the head.
5. When the weight is detected, the value flashes 3 times on the display and then it remains displayed.
6. After measurement, press and hold the "ON/OFF - TARE" button to power off the scale. If there is no weight change, the display will continue to show the value of the last measurement for one minute and then it will automatically power off.

Tare function

This function allows to weigh and exclude from the measurement the weight of any mat or towel placed on the scale to weigh a naked baby.

1. Press the "ON/OFF - TARE" button to power on the scale. The display will show «**Hello**» (3) and after a few seconds «**0.00**» (2).
2. Place the mat or towel on the scale, the weight will flash 3 times on the display and then it will remain displayed.
3. Press the "ON/OFF - TARE" button, the display will show «**Tare**» (4) and then «**0.00**» (2).
4. Measure the weight according to point 3 in the "How to use the scale" section.

Warnings

1. The display shows «**LO**» (1): low battery, replace the batteries.
2. The display shows «**O_LD**» (5): overload, weight exceeding 20 kg, remove baby from scale.
3. The display shows «**UNST**» (6): scale placed on an unstable surface, reposition the scale on a flat, stable surface.
4. The display shows «**OUTZ**» (7): load cell damaged, contact customer service.

SAFETY, CARE AND DISPOSAL

- This device must be used exclusively as described in this manual. The manufacturer shall not be liable for any damages caused by improper use.
- This device has been manufactured using sensitive components and it must be handled with care.
- Make sure to safeguard the product from:
 - water and humidity
 - extreme temperature
 - impacts and falls
 - contamination and dust
 - direct sunlight
 - heat and cold
- Do not use the device if you deem it damaged or if you notice any anomaly.
- Never open the device.
- Remove the batteries when you do not use the device for an extended period of time.
- Never leave a baby unattended during weighing operations.
- Only use the scale with weights lower than 20 kg / 44 lb. Using the scale with heavier weights may damage the device and be hazardous.
- Do not place anything on the scale when it is not in use.

How to care for the device

Clean the device using only a soft, dry cloth.

Disposal

Batteries and electronic devices must be disposed of according to local regulations and not with domestic waste.

TECHNICAL SPECIFICATIONS

Maximum weight capacity: 20 kg / 44 lb
Minimum weight capacity: 60 g / 0.13 lb
Accuracy: 10 g / 0.02 lb
Power supply: 4 AA batteries (INCLUDED)

We reserve the right to apply technical changes.

WARRANTY

This device is covered by a **2-year** warranty from the date of purchase, against manufacturing defects, valid only upon submission of the purchase receipt. Opening or tampering with the device voids the warranty. The warranty does not cover accidental damage caused by improper use or failure to abide by the instructions for use.

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Equipment Co., Ltd